

E-connect 5

Endomotor With Built-in Apex Locator





Torque Calibration For precise torque output




OLED Screen
Digital OLED Display



Ambidextrous Display Screen can be rotated as per Left or right hand use



Powerful Battery
Heavy duty 1500mAh
Lithium battery

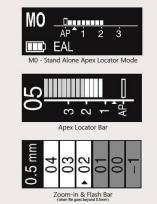

APEX LOCATOR

Apical control with rotary instruments

Position of file is displayed on screen and the motor rotates in reverse direction when file reaches pre-set apical working length, to prevent over instrumentation

Manual canal measurement

Manual file clip can be used to measure canals precisely with hand files

- Multi frequency technology with high accuracy & precision
- Apex locator flash bar with zoom-in function

Working Modes:

FWD - Clockwise Continuous Rotation

REV - Anti-clockwise Continuous Rotation

REC - Reciprocation Movement

First-in-class cordless open motor with adjustable angles of rotation ranging from 30° to 370° for reciprocation as per user's need

First-in-class feature of speed adjustment in reciprocation mode

ATC - Adaptive Torque Control

ATC manages the file stress in a smart way. When torque limit is reached, instead of going into reverse mode, motor reciprocates till the stress is reduced. Once back to below pre-set torque level, file will again rotate in continuous rotation

Integrated mode features:

- Auto Start & Stop Motor starts working when the file enters inside root canal, and stops when withdrawn
- Apical Reverse When the file reaches pre-set apical reverse point, user can define whether the file starts moving in reverse direction
- Apical Slow Down To preserve the cemento-dentinal junction and avoid file separation, speed of motor slows down when file is in apical region

Specifications:

Rotational Speed Range : 120-1000 rpm Torque Range : 0.5-4.0 N.cm Memory Programs : 11 (M0 to M10) Battery Capacity : 1500 mAh

Eghteeth