





CEFLA s.c. Via Selice Provinciale 23/a • 40026 Imola • Italy t. +39 045 8202727 • 045 583500 info@newtom.it

newtom.it

RXDC INTRAORAL X-RAY UNIT.

Outstanding quality and innovation, exceptional features.

RXDC X-VS

PERFECT.VISION

2

RXDC efficiency stems from a combination of advanced technology and an outstanding capacity to produce high definition images. The RXDC X-ray unit provides top-flight performance, practicality and technology. The RXDC features a constant potential high frequency (DC) generator and a very small focal spot (0.4 mm) capable of providing sharp, detailed images while ensuring working comfort and low doses for the patient.

Higher performance with RXDC, the X-ray unit that combines high definition imaging, ergonomic design and low X-ray doses.

О

Thanks to rectangular

collimation and the ECO

Mode parameters, the

patient exposure to

X-rays is minimal.

DOSE



SUPERIOR DIAGNOSTIC QUALITY

Obtained in just a few simple steps, all images are high resolution.



ADVANCED TECHNOLOGY

focal spot.

The NewTom RXDC high-frequency X-ray unit is based on NewTom's know-how with a 30 cm source distance and 0.4 mm



MINIMAL RADIATION **VERSATILE AND EASY TO INSTALL**

Easy, fast installation with multiple positioning options. NewTom RXDC is available in both a wall-mounted and a trolley-mounted version.





PRECISION DIAGNOSTICS.

Immediate diagnosis, excellent results.

Focal spot 0.4 mm and power 70 kV, 8 mA, highfrequency constant potential generator. Cutting-edge technology for extremely detailed images. The RXDC is extremely reliable: constant-potential design ensures image generation is unaffected by power fluctuations.



USER-FRIENDLY CONTROL

A practical, user-friendly handheld unit, designed for immediate, precise X-ray image acquisition, allows easy selection of the most suitable programme. Moreover, it allows users to control the exact emitted dose and the tube temperature via the sequential exposure graph. A wi-fi version is also available.



SUPERIOR PERFORMANCE AND TOP-CLASS ERGONOMICS.

Thanks to the protractor with graduated scale, positioning of the arms and the head is stable, effective and fully adaptable to your work. Consists of extruded aluminium arms with an integrated self-balancing system – available in the following lengths: 40 cm (15.7") – 60 cm (23.6") – 90 cm (35.4"). The adjustable wall support ensures maximum installation versatility.


Increased X-ray parallelism and an incorporated collimator allow the RXDC to achieve a source-to-skin gap of 30 cm. The RXDC provides pin-sharp, precise images with outstanding detail.

Extremely practical and versatile, RXDC can be used together with any type of direct or indirect digital sensor and X-ray film. Featuring 28 levels of sensitivity, it ensures sharp images in any situation.



The RXDC unit can also be set up with shutters and an (optional) rectangular collimator to define the body area that will be exposed and so reduce the received dose. Maximum attention to staff and patient health, while ensuring sharp, high definition image quality.



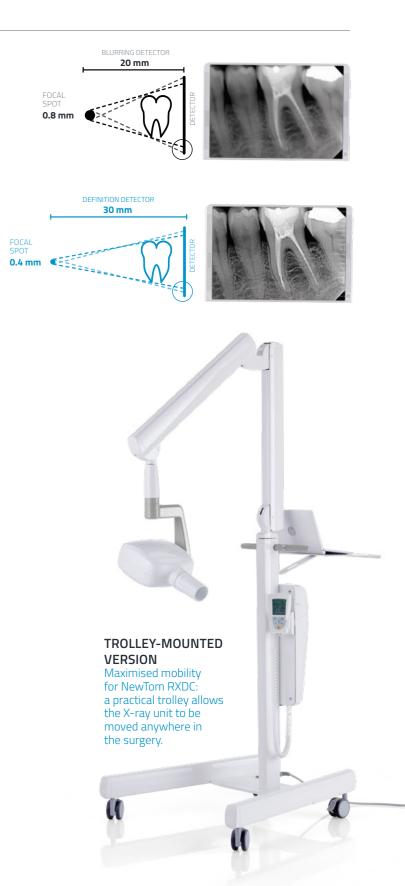

RXDC X-VS

PERFECT.VISION

X-VS INTRAORAL SENSOR.

For high quality low dose dental diagnostics.

Innovative ergonomics, direct USB plug-and-play connection, high definition and immediate results make the X-VS with HR technology the most advanced and suitable sensor for your surgery. Simplicity of use and image acquisition - combined with advanced real-time digital technology - improve quality of work.

The latest generation of X-VS image processing software aims to improve diagnostic efficacy. With excellent image resolution and an intuitive software interface, X-VS makes reading images easier and better suited to the purpose. This translates into a comfort zone personalized for each professional and for each appointment.

0

MULTI-VISION DIAGNOSTICS

High definition image acquisition in a few simple steps with filters optimised for every clinical need.



HR PLUG & PLAY DIGITAL TECHNOLOGY

FOP multilayer sensor, Rounded profiles and sturdy and reliable with ergonomic design to direct USB connection. adapt to the oral cavity. Maximised active area ensuring an extended

OPTIMAL

view.

ERGONOMICS

 \cap



INTEGRATED IMAGE MANAGEMENT

The NNT software manages, processes and shares the acquired images on PC; also visible on iPad**.



CUSTOMISED DIAGNOSTICS

Available in two sizes for maximum adaptability to the dimensions of the patient's oral cavity. Excellent working comfort and positioning, ensured by ergonomic sensors with rounded corners. A set of innovative filters allows customised tests to be carried out to improve the diagnostic vision.

WATERPROOF WITH IP67 PROTECTION RATING

RELIABLE AND ERGONOMIC.

Multivision for real-time quality diagnosis. The X-VS intraoral sensor offers extraordinary performance, practical ergonomics and high technology, offering a perfect balance between comfort and cutting-edge technology. X-VS is impact-resistant, dust-resistant, IP67 certified (water-resistant) and can be used with all X-ray systems.

X-VS means real-time diagnostics, direct USB plug-and-play connection, high definition and immediate results. X-VS uses iRYS, the all-in-one software ideal for diagnostics, communication and management of intraoral imaging: perfect for storing, managing and printing images in perfect synchronism with any other devices already in the surgery.

INNOVATIVE ERGONOMICS

Ergonomic design, rounded corners and a flexible lead make the X-VS a practical, ergonomic and intelligent sensor. This speeds up the work and makes it more practical, maximising patient comfort. Designed to adapt perfectly to the anatomy of the oral cavity, X-VS maximises both the active area and positioning comfort. Ergonomic positioners ensure optimal sensor placement.

X-VS maintains a perfect combination of first-rate comfort and cutting-edge technology Patient comfort is ensured by ergonomics and automatic acquisition, helping real-time diagnostics: it also allows the dentist/assistant to be always next to the patient for an uninterrupted workflow.

With X-VS the captured images are immediately displayed. Quick and easy sharing, communication and storage for an optimised workflow. Following acquisition, images are loaded directly onto the PC.

From here they can be consulted, printed and shared via the iPad App** or a free image viewer.

NEV The in Thank a set o Equipp MultiV visual MU Four-I colum interce The Fi


NEWTOM ADAPTIVE MULTIVISION

The innovative NewTom **ApT** (Adaptive Picture Treatment) filters have been specially developed to meet the needs of professionals. Thanks to proprietary algorithms optimized for the X-VS sensor, this function allows users to simultaneously acquire, display and share a set of images (up to 5), each with a specific improvement useful for highlighting anatomical details with different levels of sharpness. Equipped with the powerful NNT software, X-VS now allows more advanced and versatile image processing filters to be used in **Adaptive MultiVision** mode. You can select which filters to use from pre-set families or define customised ones based on individual diagnostic or visual preferences. This enhances diagnostic efficiency.

MULTILAYER SENSOR

Four-layer sensor, Caesium lodide scintillator with column-like micro-structures that preserve image quality; intercepts the X-ray beam and converts it into visible light. The Fibre Optics Plate collimates the radiation onto the sensor and protects it against X-ray penetration. The CMOS acquisition device and the electronics convert the light into a high definition digital image



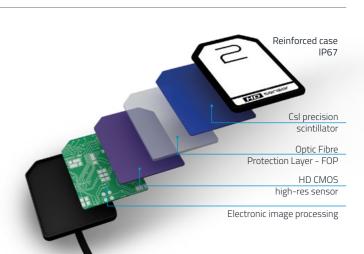

TECHNICAL SPECIFICATIONS.


RXDC X-RAY UNIT


Generator	Constant potential, microprocessor-controlled	
Working frequency	145 - 230 KHz with self-adjustment (typically 175 KHz)	
Focal spot	0.4 mm (IEC 336)	
Total filtration	2.0 mm Al@ 70kV	
Anode current	4 / 8 mA	
Voltage at X-ray tube	60 / 65 / 70 kV (*)	
Exposure times	0.020 – 1.000 seconds, R'10 and R'20 scale	
Source-skin distance	20 and 30 cm	
Irradiated field	Ø 55 mm and Ø 60 mm round	
Additional collimators	35 x 45 mm rectangular, 31 x 41 mm and 22 x 35 mm, for sensors size 2 and size 1	
Power supply	50/60 Hz, 115-120Vac ±10% or 230-240Vac ±10%	
Duty Cycle	Continuous operation with self-adjustment up to 1s/90s total	
Arms (for Standard version only)	Available in 3 lengths: 40 cm (15.7") - 60 cm (23.6") - 90 cm (35.4")	
Max. arm extension	230 cm, from wall	
Versions	Standard (wall mounted) or Mobile (on portable cart)	


(*) values depend on the country where the product is marketed


X-VS SENSOR	SIZE 1 – STANDARD	SIZE 2 – LARGE
Outside dimensions (mm)	38.9 x 24.9	41.9 x 30.4
Thickness (mm)	5.3	5.7
Pixel matrix	1500 x 1000	1700 x 1300
Pixel size (µm)	20	20
Max. resolution (lp/mm)	25	25
Grey level depth	14 bit acquisition - 16.384 maximum levels of grey	
Scintillator technology	Csl (Caesium lodide) with micro-columnar structure	
Direct exposure protection	FOP (Fibre Optic Plate)	
Degree of Case protection	IP 67 (Guaranteed against liquid or dust infiltration)	
Compatibility with X-ray generators	Any AC or DC technology X-ray generator with kV values in the 60 – 70 kV range and precision control of exposure times	
Connectivity	Direct USB to PC	
Acquisition software (for PC)	iCapture with dedicated filters for third party software	
Image management software (for PC)	NNT (compliant with ISDP©10003:2018 in accordance with EN ISO/IEC17065:2012 - certificate number 2019003109-1) and iPad NNT viewer App** (free)	
Supported protocols	DICOM 3.0, TWAIN, VDDS	
DICOM nodes	IHE compliant (Print; Storage Commitment; SR document; WorkList MPPS; Query Retrieve)	
MINIMUM SYSTEM REQUISITES		
Supported operating systems	Microsoft® Windows® 10 Professional 64 bit	
Display settings	1280 x 1024; 1344 x 768 or greater, 16 million colours	
Port	USB 2.0 or later versions	
Power supply	5 VDC, 500 mA (via USB)	

Specifications subject to change without prior notice.

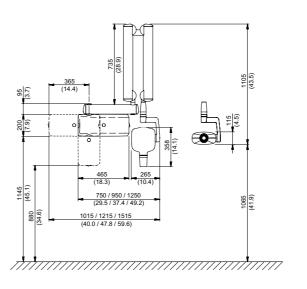

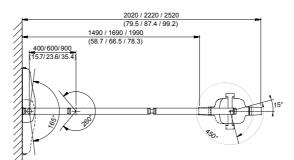

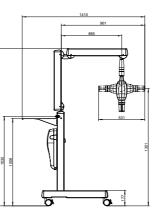


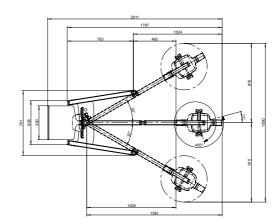

1145

Available on the App Store



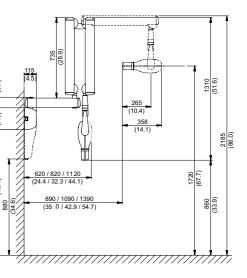

Dimensions in centimetres (dimensions in inches)



11