


CATALOG 2018

● Amalgamators

TAC UltraMix

LINEA TAC

AMALGAMATOR

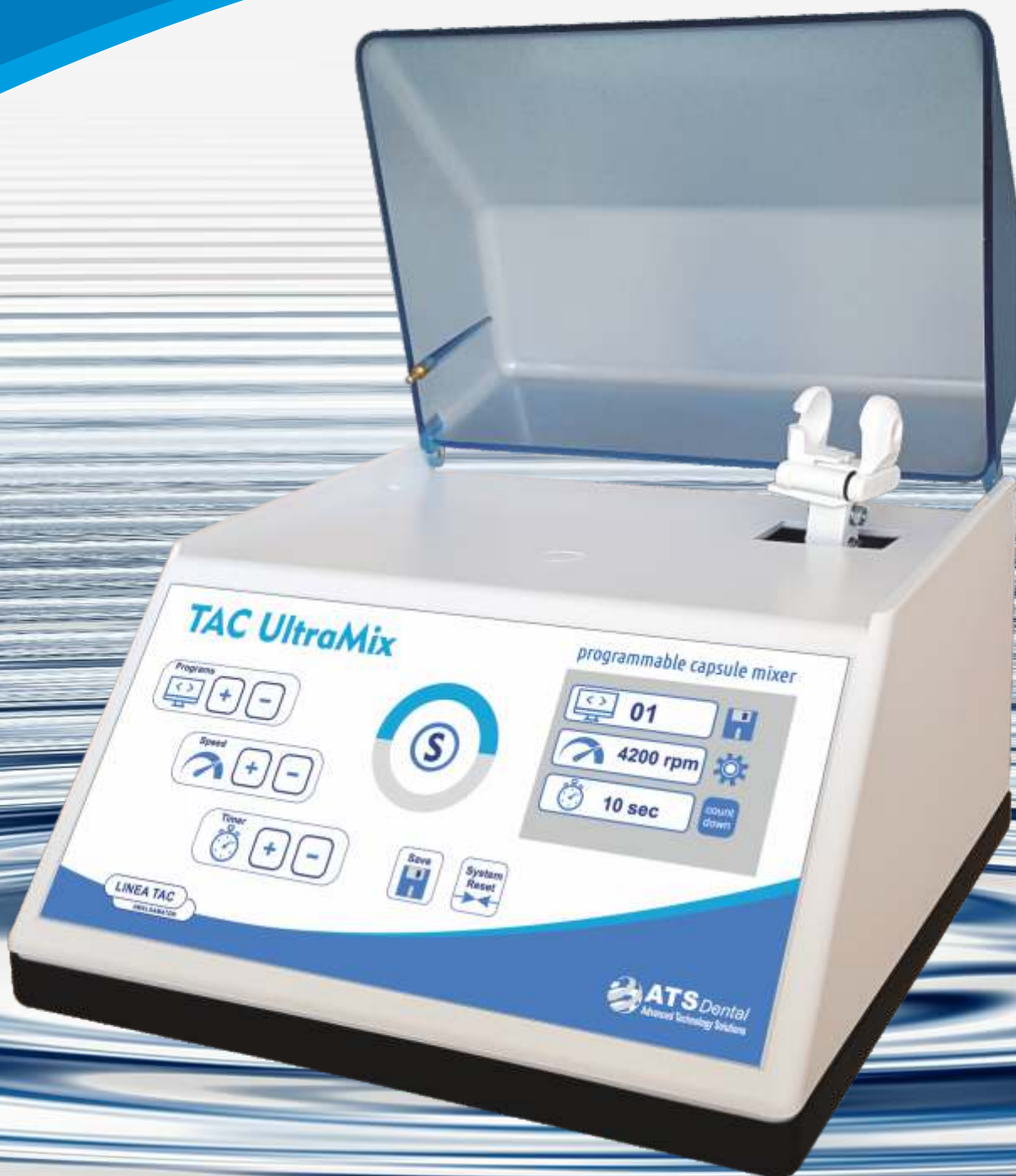

Mixing.Strong.Always

TAC UltraMix

LINEA TAC

AMALGAMATOR



fully computerized and programmable capsule mixer



TAC UltraMix is a modern universal mixing device for amalgam, glass ionomer cements and other predosed dental materials in capsules and injection capsules. This well-engineered and powerful unit has a reliable **«speed and time control»**. The microprocessor controls:

- **«10 Programs»**
- **«Speed»** that can be adjusted **from 2600 rpm to 6000 rpm**
- **«Time»** that can be adjusted **from 0 to 99 seconds**.

Cycles per minute are continuously verified.

«S»
- mixing Start

«Programs»


- Scrolls through the program
- Choice of 10 programs

Speed

- Sets the rpm. Motor speeds from 2600 to 6000 rpm

Timer

- Sets the time: from 1 to 99 sec.

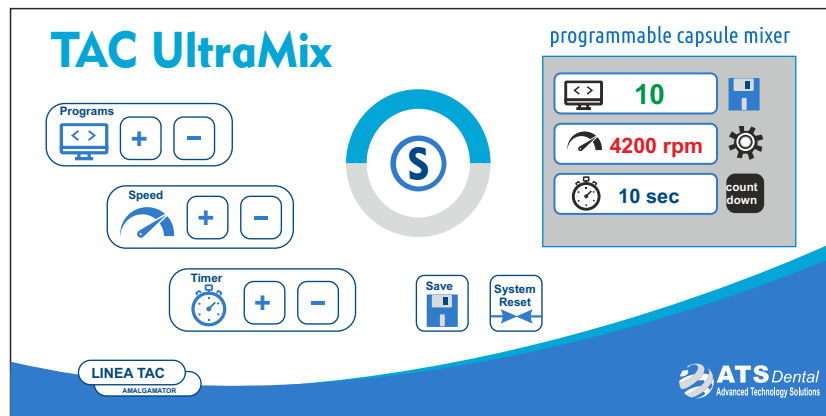

«Graphic Display»

- clearly displays program in use, speed and time

«Save»

- save values of the programs 1-10

«System Reset»

- set all the values to the default parameters

Simply to use

The touch control system is quick and simple to operate. The capsule holder design makes loading and removal easy

Touch Pad Selector

The 2600 - 6000 rpm speed range ensures the most complete mixing possible

Universal Arm Design

Multi use mixer suitable for all types of encapsulated dental materials

Reinforced Carbon-Fiber Arms

More durable than metal arms

Vibration absorbers. Quite operation and low vibration

Minimal vibration of all external parts ensures it will not move from its standing position

Protective Cover

Covers arms and capsule during high speed mixing to reduce noise. It will immediately stop mixing when the protective cover is open

Ramp to seed quickly

It goes to 0 to set cycles per minute immediately

Easy to disinfect

The hygienic one surface covers facilitate infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.

Model	Code	Description
TAC-UltraMix	TAC UTX	Fully computerized and programmable capsule mixer

TECHNICAL DATA

Model	TAC-UltraMix	Class MD	1 rule 12 in Annex IX 93/42/EEC
Dimensions	300x220x170	Console	IP21
Materials	ABS	7488 Classification	Type 2
Noise	< 65 dBA	Operating Mode	30° ON – 180° OFF
Supply Voltage	AC 100-240V out 30V	Max Working Time	30'
Frequency	50/60 Hz	Fuses	1xT2.0A 250V

Reinforced Carbon-Fiber Arms



TAC 400-M

LINEA TAC

AMALGAMATOR



Mixing.Strong.Always

TAC 400-M

LINEA TAC

AMALGAMATOR

Features and Benefits

computerized and programmable capsule mixer

TAC 400-M is a computerized and programmable capsule mixer with 10 pre-set memories. The microprocessor controls to provide extreme accuracy for a perfect mix-to-mix consistency and performance, continuously verifying cycles per minute.

- «10 Programs»
- «Speed» **4200 rpm**
- «Time» that can be adjusted **from 1 to 99 seconds**.

«LED Display»

- displays program in use and time

«M»

- changes 10 programs

«+ -»

- changes mixing times

«S»


- mixing Start

Reinforced Carbon-Fiber Arms



TAC-400M

programmable capsule mixer

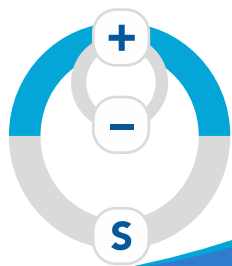

M

LINEA TAC

AMALGAMATOR

ATS Dental
Advanced Technology Solutions

Simply to use

The touch control system is quick and simple to operate. The capsule holder design makes loading and removal easy

Touch Pad Selector

4200 rpm speed range ensures the most complete mixing possible

Universal Arm Design

Multi use mixer suitable for all types of encapsulated dental materials

Reinforced Carbon-Fiber Arms

More durable than metal arms

Vibration absorbers. Quite operation and low vibration

Minimal vibration of all external parts ensures it will not move from its standing position

Protective Cover

Covers arms and capsule during high speed mixing to reduce noise. It will immediately stop mixing when the protective cover is open

Ramp to speed quickly

It goes to 0 to 4200 cycles per minute immediately

Easy to disinfect

The hygienic one surface covers facilitate infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.



Model	Code	Description
TAC 400-M	TAC 1003	Computerized and programmable capsule mixer

TECHNICAL DATA

Model	TAC-400M	Class MD	I rule 12 in Annex IX 93/42/EEC
Dimensions	300x220x170	Console	IP21
Materials	ABS	7488 Classification	Type 1
Noise	< 65 dBA	Operating Mode	30" ON – 180" OFF
Supply Voltage	AC 230V	Max Working Time	30'
Frequency	50/60 Hz	Fuses	1xT2.0A 250V



TAC 200-S

LINEA TAC

AMALGAMATOR



Mixing.Strong.Always

TAC 200-S

LINEA TAC


AMALGAMATOR

Features and Benefits

computerized capsule mixer

TAC-200S is a computerized capsule mixer. The microprocessor controls to provide extreme accuracy for a perfect mix-to-mix consistency and performance, continuously verifying cycles per minute.

- «Speed» **4200 rpm**
- «Time» that can be adjusted **from 0 to 30 seconds**.

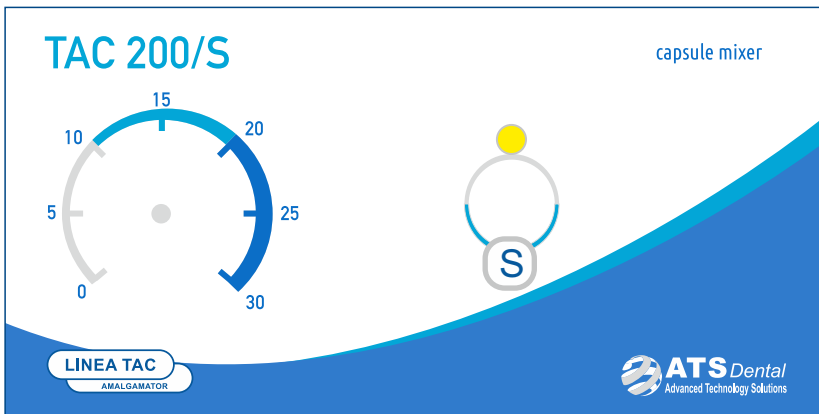

«Analogic time setting»

- It changes the time of mixing

«S»

- mixing Start

Reinforced Carbon-Fiber Arms



Simply to use

The analogic/touch control system is quick and simple to operate. The capsule holder design makes loading and removal easy

Touch Pad Starting

The 4200 rpm speed range ensures the most complete mixing possible

Universal Arm Design

Multi use mixer suitable for all types of encapsulated dental materials

Reinforced Carbon-Fiber Arms

More durable than metal arms

Vibration absorbers. Quite operation and low vibration

Minimal vibration of all external parts ensures it will not move from its standing position

Protective Cover

Covers arms and capsule during high speed mixing to reduce noise. It will immediately stop mixing when the protective cover is open

Ramp to speed quickly

It goes to 0 to 4200 cycles per minute immediately

Easy to disinfect

The hygienic one surface covers facilitate infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.



Model	Code	Description
TAC 200-S	TAC 1032	Computerized capsule mixer

TECHNICAL DATA

Model	TAC-200S	Class MD	I rule 12 in Annex IX 93/42/EEC
Dimensions	300x220x170	Console	IP21
Materials	ABS	7488 Classification	Type 1
Noise	< 65 dBA	Operating Mode	30° ON – 180° OFF
Supply Voltage	AC 230V	Max Working Time	30'
Frequency	50/60 Hz	Fuses	1xT2.0A 250V





LINEA TAC

AMALGAMATOR

Acquisition of «LINEA TAC AMALGAMATOR»

LINEA TAC

AMALGAMATOR

CELEBRATING
40
YEARS

ATS DENTAL, manufacturing motors used in endodontics and implantology, has expanded its operations with the acquisition of the «Linea Tac Amalgamator business» in March 2016. Acquiring «Linea Tac» is a significant step forward in achieving a leadership role in the field of Amalgamators. Linea Tac

Via Vecchia Provinciale Lucchese, 49/F

Zona Industriale "Mauro Maltinti"- Loc. Masotti

51030 SERRAVALLE PISTOIESE (pt)

Tel: +390573518137 - Fax: +3905731711140

www.atsdental.it - info@atsdental.it



Import Dealers Welcome !

Printed in Italy - FEB 2017

ATS DENTAL srl reserves the right to modify this catalog without prior notice

