

new|D|A INTRAORAL SENSOR


THE EVOLUTION OF YOUR IMAGE



Size 1
active area:
20x30 mm



Size 2
active area:
26x36 mm

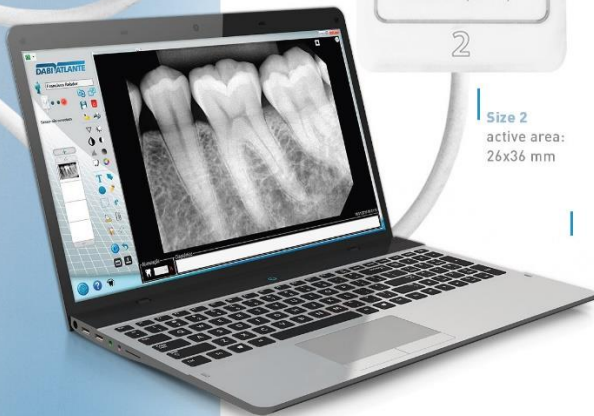

The intraoral sensor New IDA provides efficient images for a sharper, accurate and thus safer diagnosis.

The New IDA connects you to the future, optimizes your intraoral radiography workflow and manages your patients' images efficiently.


Save time, money, perform accurate diagnoses, impress your patients and help the environment!

The New IDA is compatible with most X-Ray equipments available on the market. The software is user-friendly and can also be installed on as many computers as you wish.

It's time to leave analogic behind and digitalize your clinic.

NEW IDA BENEFITS

- Easy access among colleagues
- Edit and save your images
- Environmentally friendly
- Accurate diagnoses
- Reduce film and processing costs
- Optimize your time
- Impress your patients
- Low radiation levels

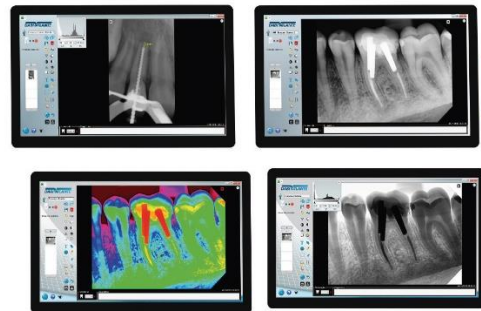

DABI ATLANTE

FEATURES

Compatibility

The intraoral sensor New IDA is compatible with most X-Ray equipments available on the market.

Environmentally friendly

Put an end to your chemical waste. Do your part and contribute to the environment.

Unlimited licenses

The New IDA comes with no license limits and you can install its software on any computer. Save your images in a network folder and access it from your device anytime you want.

Share your images and reports

Save your images in a variety of file formats, including TIFF, BITMAP, PNG, DICOM and JPEG. Create PDF reports and share them swiftly with your colleagues and patients via e-mail.

Save time

While the analog systems takes up to 5 minutes to generate images, it takes New IDA only 5 seconds.

Impress your patients

New IDA takes your clinic into the future. Your patients will be simply amazed by its high technology. Demonstrate your work to patients with a wealth of details, comparing before and after images.

Organized images

Register your patients, capture images, organize and report all your images in just a few seconds. With the New IDA you will be able to access and share any images in a timely manner.

Easy cable replacement

The New IDA uses an easy to replace 2.0 USB cable (3 meters long) in addition to a internal protective system.

Save money


Stop wasting your money with film and chemicals. New IDA is the most cost-effective investment.

Sharper images

Capture high-definition images in real time. The New IDA has a 25 line pairs per millimeter resolution and also advanced image processing applications, generating sharper images for the most accurate diagnosis.

Full-featured

The New IDA offers a wide range of user-friendly tools, such as: rotation, zoom, brightness, contrast, positive/negative, color adjustment, text insertion fields, arrows, circles, cropping, selection, measurement, histogram adjustments, odontogram and full screen image diff tool.

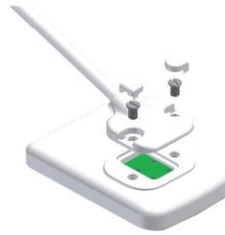

TECHNICAL INFORMATION

Operational System: Windows 7 Home, Start, Pro (x86 and x64); Windows 8 Pro - (x86 and x64), Windows 8.1 Pro - (x86 and x64); Windows 10 Home and Pro - (x86 and x64); CPU Intel Core i3 3.0 GHz or higher; RAM 4GB DDR2 SDRAM; HDD 500GB or higher; PCI - ; NIC - ; Video - ; USB Connector 2x USB 2.0; Monitor 1024 x 768; Power Supply - .

DABI ATLANTE

INTRAOURAL
SENSOR

new IDA